

# Rehab Hock Protector – Buying Guide



## KRUUSE Rehab Hock Protector

The hock protector can be used before and after surgery or in connection with chronic arthritic conditions such as osteoarthritis, OCD, ligament injuries, arthritis and other problems.

A well-adapted joint protector that repels water is important, as the joint protector comes near the ground and therefore easily gets wet.

A wet joint protector cools down the joint and impairs the healing process.


The KRUUSE Rehab Hock Protector is made of neoprene and is waterproof. The inside of the joint protector will remain dry even under extremely poor weather conditions. The hock protector is placed with the point of the hock out through the opening in the protector and then stretched around the leg proximally and distally to the joint.

Size	Cat. No	Suitable for
XXX-Small	279933	Chihuahua, Japanese Chin, Maltese, Pomeranian, Japanese Terrier, Yorkshire Terrier, Papillon
XX-Small	279932	Bichon Frisé, Coton de Tulear, Tibetan Spaniel, English Toy Spaniel, Yorkshire terrier
X-Small	279839	Sheltie, Cocker Spaniel, Fox Terrier, French Bulldogs, Yorkshire terrier, Miniature Schnauzer, Cavalier king Charles Spaniel, Shih Tzu, Pomeranian
Small	279840	Kelpie, Whippet, Bulldog, Beagle, Dachshund, Pembroke Welsh Corgi, Boston Terrier
Medium	279841	Labrador, Border Collie, Springer Spaniel, Labrador Retriever, Golden Retriever, Poodle, Boxer, Pointer (German Shorthaired), Doberman, Australian Shepherd
Large	279842	German Shepard, Boxer, Rhodesian Ridgeback, Siberian Husky, Belgian Malinois, Alaskan Malamute
X-Large	279884	Bernese Mountain Dog, Great Dane, Leonberger, Irish Wolfhound, Mastiff, Saint Bernard

### User Manual

The hock protector is placed with the tip of the hock out through the opening in the protector and then stretched around the leg above and under the joint.

# Rehab Hock Protector – Size Guide

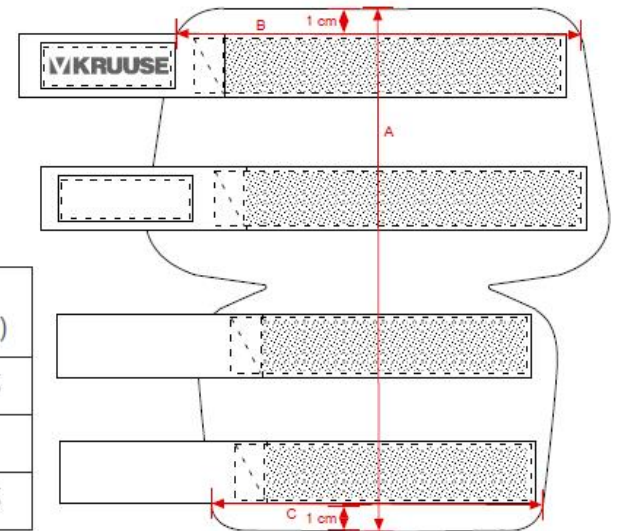

	Size chart KRUUSE Rehab Hock Protectors in cm	XXXS (cm)	XXS (cm)	XS (cm)	S (cm)	M (cm)	L (cm)	XL (cm)
A	Total height, from edge to edge	8	12	14	16	17.5	19	20.5
B	Total upper width 1 cm from edge	6	8	10	12	14	15	16
C	Total lower width 1 cm from edge	4	5.5	7.5	9.5	11.5	12.5	13.5

# Rehab Knee Protector – Buying Guide



## KRUUSE Rehab Knee Protector

Designed by ReDog of Sweden™

*"The knee joint is an advanced joint, which often gives the dog long-term convalescence with the risk of subsequent problems following the injury."*

*When the dog's knee joint becomes injured, the musculature quickly weakens and the joint becomes unstable, which can lead to meniscus injuries and arthrosis changes in the articular cartilage. With the KRUUSE Rehab Knee Protector you can help the dog directly in conjunction with the diagnosis being made and support the joint and musculature until it is no longer deemed as being necessary.*

*My experience is that dogs with knee joint protection before and after the operation of the front cruciate ligament always return to normal function more quickly after injury.*

*The Knee Protector is constructed so that it gives as equal stabilisation as possible in order to imitate the front cruciate ligament's normal function, which makes the dog feel safe in using the protector and often it goes better with the protection in just 5 minutes compared with what they did before the diagnosis was made."*

*Stefan Rosén, ReDog AB*

The knee protector can be used before and after surgery and in cases of chronic arthritis, cruciate ligament injuries, arthritis, OCD, meniscus injury, kneecap problems, tendinitis etc.

## User Manual


The knee protector must be pulled up along the leg to sit well. This is done by stretching the knee protector up on the dog's leg to cover the dog's knee properly. Lead the soft strap that holds the protector in place under the abdomen and over the back and fasten it on the corresponding touch fastener on the outer side of the protector.

## Care

Wash by hand and air-dry.

# Rehab Knee Protector – Size Guide

Size	Cat. No	Suitable for
XXX-Small, left	279944	Chihuahua, Japanese Chin, Maltese, Pomeranian, Japanese Terrier, Yorkshire Terrier, Papillon
XXX-Small, right	279945	
XX-Small, left	279942	Bichon Frisé, Coton de Tulear, Tibetan Spaniel, English Toy Spaniel, Yorkshire terrier
XX-Small, right	279943	
X-Small, left	279849	Sheltie, Cocker Spaniel, Fox Terrier, French Bulldogs, Miniature Schnauzer, Cavalier king Charles Spaniel
X-Small, right	279850	
Small, left	279851	Kelpie, Whippet, Bulldog, Beagle, Dachshund, Pembroke Welsh Corgi, Boston Terrier
Small, right	279852	
Medium, left	279853	Labrador, Border Collie, Springer Spaniel, Labrador Retriever, Golden Retriever, Poodle, Boxer, Pointer (German Shorthaired), Doberman, Australian Shepherd
Medium, right	279854	
Large, left	279855	German Shepard, Boxer, Rhodesian Ridgeback, Siberian Husky, Belgian Malinois, Alaskan Malamute
Large, right	279856	
X-Large, left	279857	Rottweiler, Bullmastiff, Bloodhound, Dogue de Bordeaux, Cane Corso, Akita
X-Large, right	279858	
XX-Large, left	279885	Bernese Mountain Dog, Great Dane, Leonberger, Irish Wolfhound, Mastiff, Saint Bernard
XX-Large, right	279886	

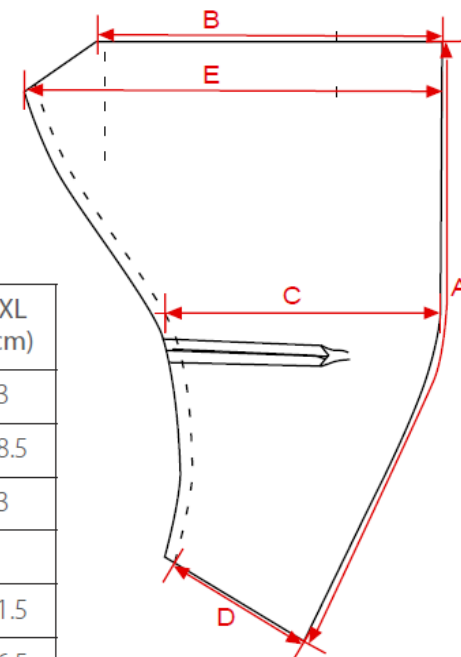

	Size chart KRUUSE Rehab Knee Brace in cm	XXXS (cm)	XXS (cm)	XS (cm)	S (cm)	M (cm)	L (cm)	XL (cm)	XXL (cm)
A	Total length, along fold edge	7	12	20.5	25	26.5	26.5	33	33
B	Width, top edge	6.5	7.5	11.5	13.5	14.5	15.5	17.5	18.5
C	Width, (0.5 cm)* 1 cm above knee dart	5*	6*	8	10	10.5	11.5	12	13
D	Width, bottom hem	4.5	5	5	5.5	6	6.5	7	8
E	Width, widest point	7	9	15	16.5	17.5	18.5	20.5	21.5
	Strap full length	50	55	60	61	90	91	96	96.5



# Rehab Elbow Protector – Buying Guide



## KRUUSE Rehab Elbow Protector

Designed by ReDog of Sweden™

*"Unfortunately, many dogs are afflicted by articular cartilage injuries or of articular cartilage wear on their elbows. Often it is related to the limited ability for the dog to be able to cope with exercise, training or competition.*

*When the dog appears to be troubled by its elbow joints, the most common symptom is lameness or stiffness, which they can then ease by moving about a bit. When dogs exhibit these symptoms, an elbow protector is the ideal solution in connection with the end of a training session in order to exert counter-pressure on subsequent swelling. Furthermore, the protection also contributes with increased heat, thus providing the muscles with better oxygen supply and better circulation conditions.*

*For many dogs, the elbow protector can be decisive in how well they can adapt to their joint problems. It counteracts pain, swelling and gives muscles a chance to rest correctly."*

*Stefan Rosén, ReDog AB*

The Elbow Protector can be used before and after surgery and in connection with chronic arthritis i.e. Osteoarthritis, OCD, ligament injuries, arthritis, ED and other assessments.

### User Manual

The elbow protector is fastened above and under the joint with the aid of two touch fastener bands that run through loops in the protector so that it does not slide down the leg.


### Care

Wash by hand and air-dry.

# Rehab Elbow Protector – Size Guide



Size	Cat. No	Suitable for
XXX-Small	279931	Chihuahua, Japanese Chin, Maltese, Pomeranian, Japanese Terrier, Yorkshire Terrier, Papillon
XX-Small	279930	Bichon Frisé, Coton de Tulear, Tibetan Spaniel, English Toy Spaniel, Yorkshire terrier
X-Small	279834	Sheltie, Cocker Spaniel, Fox Terrier, French Bulldogs, Yorkshire terrier, Miniature Schnauzer, Cavalier king Charles Spaniel, Shih Tzu, Pomeranian
Small	279835	Kelpie, Whippet, Bulldog, Beagle, Dachshund, Pembroke Welsh Corgi, Boston Terrier
Medium	279836	Labrador, Border Collie, Springer Spaniel, Labrador Retriever, Golden Retriever, Poodle, Boxer, Pointer (German Shorthaired), Doberman, Australian Shepherd
Large	279837	German Shepard, Boxer, Rhodesian Ridgeback, Siberian Husky, Belgian Malinois, Alaskan Malamute
X-Large	279883	Bernese Mountain Dog, Great Dane, Leonberger, Irish Wolfhound, Mastiff, Saint Bernard

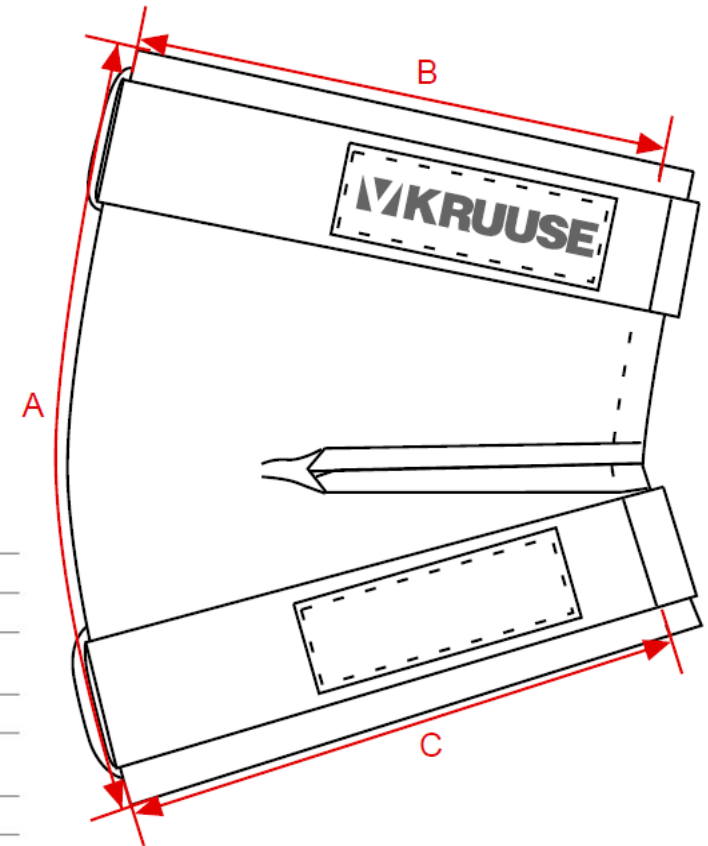

	Size chart KRUUSE Rehab Elbow Protectors in cm	XXXS (cm)	XXS (cm)	XS (cm)	S (cm)	M (cm)	L (cm)	XL (cm)
A	Total height, along fold edge	6.5	7.5	8.5	9.5	11.5	13.5	15.5
B	½ upper width, from fold edge to stitch	5.5	6.5	7.5	8.5	9.5	11	12.5
C	½ lower width, from fold edge to stitch	5	6	7	7.5	8.5	10	11.5